
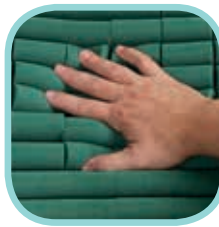


PressureGuard[®] Span-Care[®] Convertible

made with
LifeSPAN[™]

bleach-resistant stretch fabric

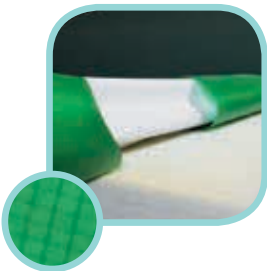

Performance, comfort and flexibility — with or without power!

The Span-Care[®] convertible air therapy surface is engineered to empower busy caregivers with an ideal blend of high-end performance and trouble-free ease of use.

In non-powered mode, use the Span-Care convertible for pressure injury prevention in at-risk patients, early intervention of existing pressure injuries, or when a step-down surface is called for. Its self-adjusting system of synchronized air cylinders and pressure relief valves provide the correct level of immersion and envelopment for any user.

Add the simple, digital control unit when powered treatment modalities are necessary. Adjustable comfort settings, audible low pressure alarm and a timed auto-firm mode ensure maximum usability.

In either mode, the Span-Care convertible delivers an unmatched combination of extreme patient comfort, crucial safety and stability, proven pressure management and maintenance-free reliability.

¹In situations where positioning on the pressure injury cannot be avoided (e.g. the individual has multiple injuries on multiple surfaces), limit the amount of time the individual is positioned on the injury.
Source: NPAUP/EPUAP Pressure Ulcer Prevention and Treatment, Clinical Practice Guideline, p 66-67

Therapy Provided

- ✓ Non-powered (Air) Reactive Pressure Redistribution
- With control unit added:
- ✓ Alternating Pressure
- ✓ Powered Pressure Redistribution

QUICK FACTS

Indications:

- Treatment of Stage 1-4 pressure injuries
- Patient can't/won't reposition frequently
- Can be positioned off the pressure injury in at least 2 positions¹
- Prevention for high risk patients

Key Advantages:

- Minimized shearing through Geo-Matt[®] surface geometry and Shear Transfer Zones[®] on bi-directional stretch cover
- Minimized risk of falls and bed entrapment through Safety Edge[™] bolster design
- Flexible treatment modalities, including alternating pressure and powered flotation
- Easy-to-use, easy-to-connect control unit

Item #	Description
SC753629	75"L x 36"W x 7"H
SC803629	80"L x 36"W x 7"H
SC843629	84"L x 36"W x 7"H
SC85V29	Span-Care Convertible for recessed deck frames
SC804229	80"L x 42"W x 7"H
SC844229	84"L x 42"W x 7"H
Control Unit: 5900	Add-on control unit
Weight limit:	Supports up to 500 lbs. (400 in powered mode) 42" models: 700 lbs. (600 in powered mode)
Weight:	Mattress: approximately 25 lbs. Control unit: 6.3 lbs.
Warranty:	Mattress: 5 years, non-prorated Control unit: 18 months, non-prorated
UL listed:	Yes
Flammability:	All models conform to NFPA 101 (Life Safety Code); CalTech bulletins # 117, 129, and 603; CPSC 16 CFR Parts 1632 and 1633

