



# Lympha-Flow

LYMPHEDEMA THERAPY

## LF900 - Superior Features and Quality - Excellent Value!

### LF900 Pump Specifications (E0651) 4 Chamber Intermittent Compression

- Adjustable Pressure Range: 10-180 mmHg (10 mmHg increments)
- **Pressure Setting Lock/Unlock function** (per Rx & safety)
- **Illuminated Start/Pause/Stop Button**
- Six Interval Speed Options
- **Therapy Timer**
- **Total Treatment Time Counter w/ Reset**
- **Continuous Treatment Mode** (for extended therapy at lower settings)
- **Four Treatment Mode Options (A,B,C,D-see diagram)**
- **Chamber Control Options: On/FOCUS, Off/SKIP (Foot, Calf, Knee & Thigh)**
- **Digital Interface, Controller & Pressure Sensors**
- Device Size: 260(W) x 180(D) x 130(H)mm (10" x 7" x 5")
- Device Weight: 2 kg (4.4 lb)
- DME Warranty: Pump-2 Years, Sleeve-6 Months

### LF900 Garments/Sleeves

(see reverse for additional info)

- **Legs (E0667)**
- **Arms (E0668)**
- **Waist**
- **Center Body**
- **Stretching Mat**
- **Convenient One-Touch Connectors**


## The LF900 System is User Friendly!

Manufactured by **DS MAREF**  
DAESUNG MAREF CO.,LTD.

## LF900


### A MODE (Pre-therapy mode)

It is mainly used to relieve the hardened inguinal lymph nodes in patients with chronic lymphedema.

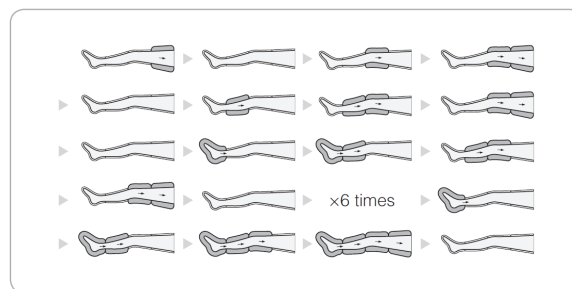


### B MODE (Lymph Mode)

This is basic mode being applied to patients with lymphedema.

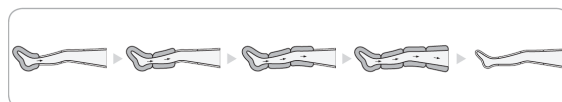


### C MODE (Distal Intensive Mode)

This is a lymphatic mode that focuses on foot area treatment for lymphedema patients with severe foot edema.

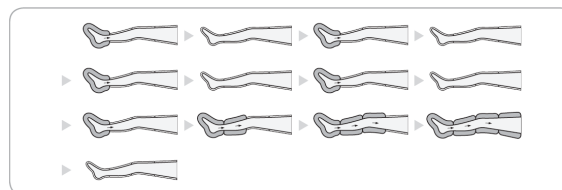

### D MODE (Soft Lymph Mode)

Patients who feel pain or uncomfortable in B-mode treatment and weakened skin and vascular tissue due to aging or chronic edema can use D-mode for gentle lymph treatment.

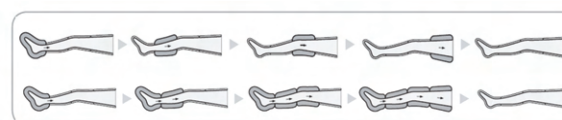

# Lympha-Flow

LYMPHEDEMA THERAPY

## LF900



### LF900 Sleeves & Garments are Specifically Designed for the LF900 System:

DS Maref has manufactured top-quality Pneumatic Compression Therapy devices since its founding in 1986. They distribute throughout the world, to 85 countries, in the medical, sports rehab, and aesthetic/spa markets with systems for both Clinical and Home use. They are a leading innovator in the design and operation of digital compression systems for treating lymphedema, phlebolymphedema, vascular disorders and wounds, DVT prophylactics, and OR tourniquets. The LF900 device is FDA approved. PDAC: E0651 (pump), E0667 (leg sleeves), E0668 (arm sleeve).

### Sleeve & Accessory Sizing Chart

#### Lower Extremity - Overlapping

SKU	Type	Size	Length	Circumference					
			Groin to Heel (Inseam)	Max Width	Mid Thigh	Mid Knee	Mid Calf	Ankle	Arch (Instep)
S401	Full Leg	Small	58cm (23")	75cm (29.5")	67cm (26.5")	65cm (25.5")	60cm (23.5")	55cm (21.5")	46cm (18")
SW401	Full Leg	Small Wide	61cm (24")	88cm (34.5")	77cm (30.5")	72cm (28.5")	62cm (24.5")	60cm (23.6")	57cm (22.5")
SXW401	Full Leg	Small Extra Wide	56cm (22")	100cm (39.5")	93cm (36.5")	85cm (33.5")	80cm (31.5")	72cm (28.5")	55cm (21.5")
M401	Full Leg	Medium Leg	71cm (28")	80cm (31.5")	72cm (28.5")	67cm (26.5")	62cm (24.5")	55cm (21.5")	39cm (15.5")
MW401	Full Leg	Medium Wide	71cm (28")	90cm (35.5")	83cm (32.5")	75cm (29.5")	70cm (27.5")	62cm (24.5")	55cm (21.5")
MXW401	Full Leg	Medium Extra Wide	71cm (28")	100cm (39.5")	95cm (37.5")	89cm (35")	84cm (33")	74cm (29")	60cm (23.5")
L401	Full Leg	Large	84cm (33")	80cm (31.5")	75cm (29.5")	69cm (27")	62cm (24.5")	56cm (22")	44cm (17.5")
LW401	Full Leg	Large Wide	86cm (34")	102cm (40")	89cm (35")	75cm (29.5")	65cm (25.5")	55cm (21.5")	44cm (17.5")
M403	Half Leg	Medium	53cm (21")			75cm (29.5")	70cm (27.5")	55cm (21.5")	46cm (18")
MW403	Half Leg	Medium Wide	53cm (21")			80cm (31.5")	75cm (29.5")	60cm (23.5")	55cm (21.5")

#### Upper Extremity - Overlapping

SKU	Type	Size	Length	Circumference					
			Shoulder to Fingertips	Chest	Axilla	Bicep	Elbow	Forearm	Wrist
M402	Arm	Medium	71cm (28")		62cm (24.25")	51cm (20")	48cm (19")	46cm (18")	32cm (12.5")
L402	Arm	Large	83cm (32.5")		72cm (28.5")	64cm (25")	56cm (22")	56cm (22")	41cm (16")
4VDRLO0	Arm/Chest	Universal	66cm (26") 84cm (33")	152cm (60") 147cm (58")	55cm (21.5")	53cm (21")	52cm (20.5")	51cm (20")	50cm (19.5")

#### Lower Extremity - Non-Overlapping

SKU	Type	Size	Length	Circumference					
			Groin to Heel (Inseam)	Upper Thigh	Mid Thigh	Mid Knee	Mid Calf	Ankle	Arch (Instep)
4LAGL5A	Full Leg	Large Non-O/L	69cm (27")	61cm (24")	55cm (21.5")	50cm (19.5")	44cm (17.5")	37cm (14.5")	33cm (13")
4LAGX40	Full Leg	XL Non-O/L	71cm (28")	75cm (29.5")	65cm (25.5")	60cm (23.5")	48cm (19")	38cm (15")	36cm (14")
4LAGY00	Full Leg	XXL Non-O/L	76cm (30")	77cm (30.5")	67cm (26.5")	60cm (23.5")	50cm (19.5")	44cm (17.5")	33cm (13")

\*Leg Extension Zippers Available - 8cm (3")

#### Upper Extremity - Non-Overlapping

SKU	Type	Size	Length	Circumference					
			Shoulder to Fingertips	Axilla to Fingertips	Axilla	Bicep	Elbow	Forearm	Wrist
G411	Arm	Universal	86cm (34")	71cm (28")	52cm (20.5")	43cm (17")			36cm (14")

#### Center Body and Waist

SKU	Type	Size	Length	Circumference				
				Upper Abdomen	Lower Thigh	Glutes	Hips	Width
4SAGF5A	Center Body	Universal	75cm (29.6")	112cm (44")	74cm (29")	152cm (60")	137cm (54")	
4ECGF01	Center Body Ext. Zippers	Universal	75cm (29.6")					10cm (4")
4WAGF1A	Waist		76cm (30")	122cm (48")				

#### Stretching Mat

SKU	Type	Chambers	Closed Width	Closed Height	Open Width	Open Length
4MDRF00	Stretching Mat	4 Chambers	61cm (24")	50cm (19.5")	61cm (24")	99cm (39")

