

MedGyn AL-106 BWL Digital Video Colposcope

Simplicity

Different from traditional optical colposcopes, the AL-106 BWL is a digital colposcope that is battery operated with wireless functionality. The colposcope has an LED light source to ensure ultra-bright visualization and a long life cycle. Its unique optics provide multiple filters to ensure high quality images and improved visualization of blood vessels.



AL-106 BWL (Front view)



AL-106 BWL (Back view)

Features & Benefits

- Wireless signal transmission
- Backup Battery
- 440K Pixels and 54X Magnification
- Multi-bright white LED lights and multiple green filters
- CCD Camera and AV Connection
- Auto/Manual focus and Image freeze function
- Digital Video Colposcope Imaging System Software



AL-106 BWL Rolling Base




Digital Video Colposcope Imaging System Software

Designed for use with MedGyn AL-106 BWL Digital Video Colposcope, the Digital Video Colposcope Imaging System Software allows for the capture and management of colposcopy images into a self-contained patient record on an easy, user-friendly and secure system.

Features & Benefits

- Complete system security, internal password protection
- Easily customizable to specific clinical requirements and preferences
- Includes a large offering of standard examination and diagnosis templates
- Full image and video capture capabilities
- Database management system
- Image measurement & annotation
- System automatically creates patient report
- Software includes: Wireless Receiver, USB Security Drive, and Footswitch

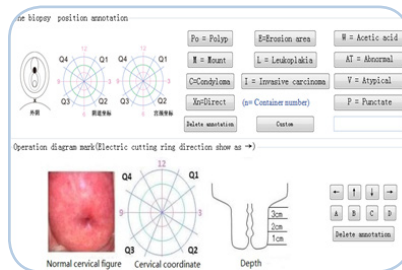

AL-106 BWL Colposcope

010471 AL-106 BWL Digital Video Colposcope with Software

010470 AL-106 BWL Digital Video Colposcope with Software and Laptop

TRUSTED PRODUCTS FOR PREMIER WOMEN'S HEALTHCARE

MEDGYN PRODUCTS, INC.

100 W. Industrial Rd. Addison, IL 60101 USA

t: +1 630.627.4105 toll-free: 800.451.9667

f: +1 630.627.0127

e: info@medgyn.com

w: medgyn.com

AL106BWL-PF1

