


# Protekt® O2 Cushion

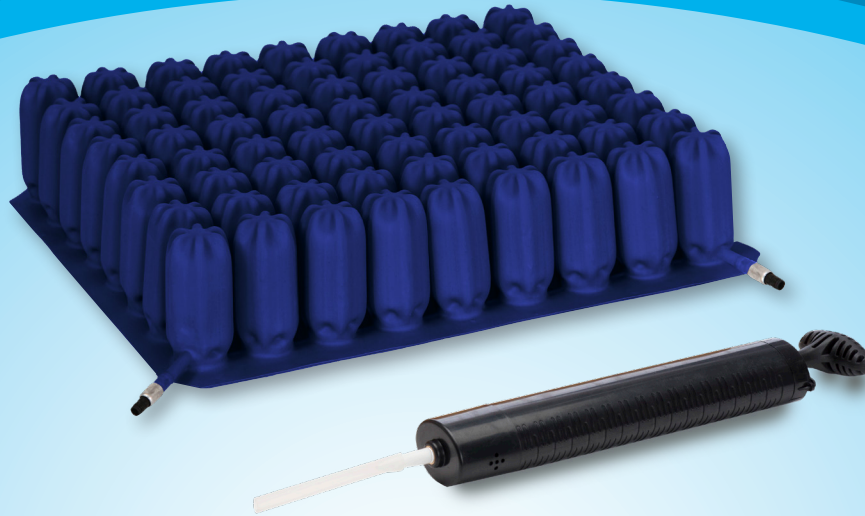

HCPCS  
E2624/  
E2625

## FEATURES AND BENEFITS:

The Protekt® O2 wheelchair cushion is an excellent choice for the prevention and treatment of pressure sores. The soft interconnected air cells allow for deep immersion providing constant low pressure across the entire seating surface.

- Adjustable 4" or 2" high air cells offers deep immersion to provide excellent pressure redistribution for tissue protection and treatment.
- Dual valve compartments address seating challenges such as pelvic obliquities, sacral seating and enhances posture control.
- Both compartments can be adjusted independently to increase positioning and stability for side to side control.
- Nylon breathable cover with non-skid bottom and safety straps.
- Easy to use hand held inflation pump.
- Flame resistant, latex free neoprene.
- Patch kit included.
- Weight Capacity: 500 lbs.
- Warranty: 5 Year Non-Prorated.

78001	16"x16"x2"
78002	18"x16"x2"
78007	18"x18"x2"
78003	20"x16"x2"
78011	20"x18"x2"
78012	22"x18"x2"
78013	22"x20"x2"

78004	16"x16"x4"
78005	18"x16"x4"
78008	18"x18"x4"
78006	20"x16"x4"
78009	20"x18"x4"
78014	22"x18"x4"
78015	22"x20"x4"



270 Washington Street, Mount Vernon, NY 10553

Tel: 855-BE-PROACTIVE (855-237-7622) • Fax: 845-205-6005

www.proactivemedical.com • sales@proactivemedical.com