

Width
74cm



160cm


Paediatric

Width
74cm



195cm

Adult

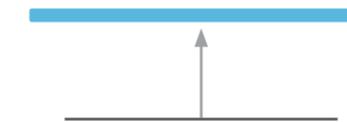

Height range 41cm to 94cm

Shower Trolleys



Designed to provide transport to and from the showering area, these mobile Shower Trolleys feature a foam padded PVC liner for added patient comfort, a 125mm wheel system with central locking & steering facility, as well as fold-down side rails for patient transfer. Paediatric models have an overall length of 160cm, and both Paediatric and Adult models are available in either hydraulic or electric height elevation. Electric models are supplied with 2 rechargeable batteries.

Models

ST7700	Hydraulic, Adult
ST7800	Electric Battery Operation, Adult
ST7600	Hydraulic, Paediatric
ST7500	Electric Battery Operation, Paediatric



Charger unit, provided with ST7800 & ST7500

- ### Standard Features
- Safe working load and lifting capacity of 240Kg (528lbs)
 - Large wheel design, with central locking and steering facility
 - Hydraulic height adjustment via foot pedal
 - Electric battery adjustment via hand switch
 - Battery operated models include a spare battery pack and charging unit
 - Foam padded PVC liner with side supports and end panels for comfort and security
 - Supplied with foam padded PVC patient support cushion
 - Side support rails fold down to assist patient transfer
 - Low gradient tilt towards foot end for water drainage via integral waste outlet
 - Plug and drain hose supplied complete with storage bracket
 - Bumpers fitted to all corners, for protection during transportation
 - Bumpers can be removed to reduce width to access doorways
 - Lowers to wheelchair height for ease of patient transfer
 - Purpose made precision bearings on all pivot points
 - Clearance beneath frame for mobile hoist
 - Choice of Adult or Paediatric models



Model ST7700
Hydraulic Adult Shower Trolley



Model ST7800
Electric Adult Shower Trolley



Model ST7700
Shown at Minimum Height
for patient transfer, with
all side rails down.