

# Chair-Air™ Series

*Alternating Pressure Chair Cushion Systems*

## **CHAIR-AIR™ 9700GR - GERIATRIC RECLINER MATTRESS AND PUMP SYSTEM**

This unique alternating pressure mattress system allows for patients who are at high risk for pressure sores to be placed in a variety of mobile devices. Patients are no longer confined to bed for the whole day due to risks of sitting without proper pressure redistribution. Chair-Air™ systems can be integrated into a geriatric chair, recliner or standard high-back wheelchair. Twenty bladders gently alternate, providing low interface pressures in addition to comfortably supporting the patient. The Chair-Air™ 9700GR system includes the 9701 alternating pump. Chair-Air™ meets California Technical Bulletin #117 for fire retardancy. An optional fire barrier is available to meet the most stringent fire safety requirements, when requested.



Recliner not included

## **CHAIR-AIR™ 9700C - ALTERNATING AIR WHEELCHAIR CUSHION AND PUMP SYSTEM**


Six individual bladders alternate to provide the patient a superior pressure redistribution seating surface. The zippered non-skid, stretch-knit cover is anti-microbial and fire-retardant. The Chair-Air™ 9700C includes the Chair-Air™ 9701 alternating pump.



Shown in cover

## **CHAIR-AIR™ 9700CAF - ALTERNATING AIR WHEELCHAIR CUSHION AND PUMP SYSTEM**

Four individual bladders alternate to provide the patient a superior pressure redistribution seating surface while the unique molded design, with slight abductor pommel, stabilizes the hips and keeps the patient in proper mid-line alignment. The zippered non-skid, stretch-knit cover is anti-microbial and fire-retardant. The Chair-Air™ 9700CAF includes the Chair-Air™ 9701 alternating pump.

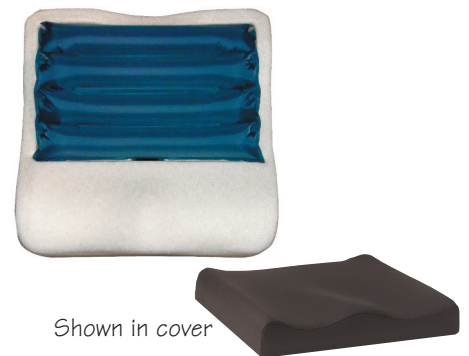

Shown in cover

## **CHAIR-AIR™ AC/DC 9701 - ALTERNATING PUMP**

The portable Chair-Air™ AC/DC Pump allows your patient to be totally mobile. A rechargeable battery backup keeps this system working up to 6 hours per charge, based on the patient's weight.



For use with Chair-Air™, Geriatric Recliner Mattress & Chair-Air™ Wheelchair Cushion

 **Blue Chip**  
Medical Products, Inc.

*"Setting the Standard"*



# CHAIR-AIR™

## SPECIFICATIONS

### CHAIR-AIR™ GERIATRIC RECLINER MATTRESS - MODEL 9702GR

**Dimensions:** 19" x 68" x 2.5 - 3.5"  
**Weight:** 6 lbs.  
**Air Cells:** Multiple - 2" - 3.5" sacs

Recliner not included


#### Mattress Cover:

Low shear, zippered multi-stretch Vyvex™ III, moisture / vapor permeable top cover, non-skid bottom. Cover is anti-microbial and meets California Technical Bulletin #117 for fire retardancy.



### CHAIR-AIR™ WHEELCHAIR CUSHION - MODEL 9702C

**Dimensions:** 18" x 16" x 3.25"  
**Weight:** 4 lbs.  
**Air Cells:** 6 - 3.25" sacs

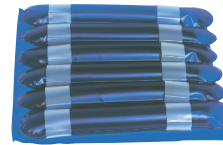

### CHAIR-AIR™ WHEELCHAIR CUSHION - MODEL 9702CAF

**Dimensions:** 18" x 16" x 3"  
**Weight:** 3 lbs.  
**Air Cells:** 4 - 2.5" sacs



#### WHEELCHAIR CUSHION COVER:

Zippered low shear stretch-knit with non-skid bottom. Cover is anti-microbial and meets California Technical Bulletin #117 for fire retardancy.

### PUMP FOR USE WITH GERIATRIC RECLINER & STANDARD WHEELCHAIR CUSHION

**Dimensions:** 14.5L" x 5W" x 12.5H"  
**Weight:** 16 lbs.  
**Electrical:** A/C 117V/60 hz D/C 12V 7 amp/hr  
**Listings:** UL CE  
**Cycle Time:** 10, 15 or 20 Minutes  
**Air Flow:** 4 - 6 Liters  
**Ground Resistance:** 0.1 ohms maximum



CHAIR-AIR™ AC/DC 9701  
ALTERNATING PUMP

#### WARRANTY

**Pumps:** 1 year limited  
**Cushion:** 1 year limited  
**Cover:** 6 months limited



7-11 Suffern Place, Suffern NY 10901

Tel: (845) 369-7535 • Fax: (845) 369-7633 • Toll Free: (800) 795-6115

www.bluechipmedical.com

"Setting the Standard"

©Blue Chip Medical Products, Inc. 2015

Alternating Pressure Mattresses  
Alternating Pressure Pads & Pumps  
Low Air-Loss Mattress Systems  
Lateral Rotation Mattress Systems  
Gel Mattresses  
Gel Mattress Overlays  
Static-Air Mattress Overlays  
Water Mattress Overlays  
Wheelchair Cushions  
Patient Lift Slings  
Side Rail Bumpers  
Cart Covers  
OR Replacement Table Pads  
Replacement Stretcher Pads  
Cubicle Curtains  
Hospital Pillows